

# Thermomac

## 6-Station Dissolution Tester

### DT Series

The measurement result is shown on the LED display. Preset and actual values can be displayed alternately.


Temperature can be set from 20 to 45 °C.

Speed can be set from 20 to 200 RPM.

Timer can be set from 1 to 999 min.

Button and Control Panel

Units and shafts are made of stainless steel.

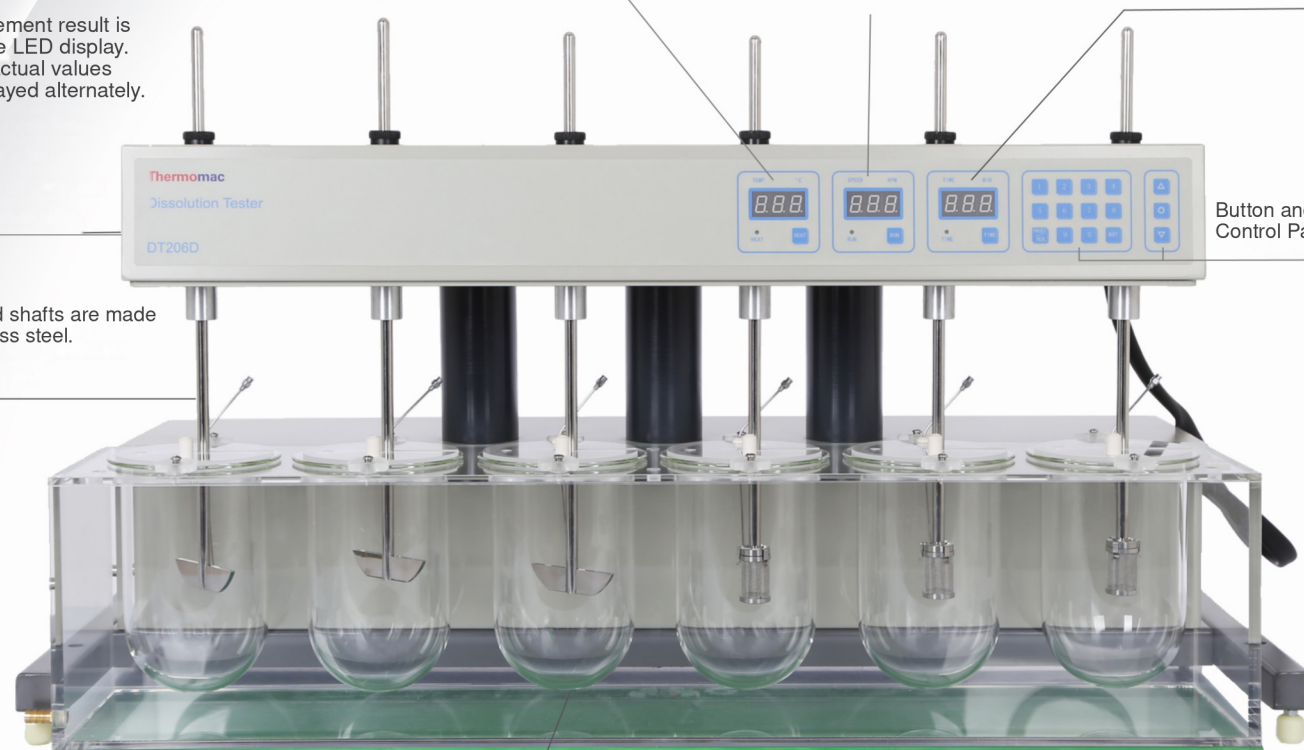

Six vessels allow samples to be tested quickly and automatically.

The magnetically pumped circular water-flow system heats evenly and maintains a uniform bath temperature.

## Technical Specifications

| Model                       | DT201                 | DT206D                 |
|-----------------------------|-----------------------|------------------------|
| Speed Range                 | 20-200 rpm            | 20-200 rpm             |
| Speed Accuracy              | ±2rpm                 | ±2rpm                  |
| Temperature Stability       | ±0.3°C                | ±0.3°C                 |
| Temperature Range (Ambient) | +45.0°C               | +45.0°C                |
| Nine-Point Preset Timing    | 5-900 min             | 1-999 min              |
| Time Accuracy               | ±0.3%                 | ±0.3%                  |
| Paddle Radial Runout        | ±0.5mm                | ±0.5mm                 |
| Basket Radial Runout        | ±1.0mm                | ±1.0mm                 |
| Power                       | 220-110V/50-60Hz/600W | 220-110V/50-60Hz/1200W |
| Dimensions                  | 240x80 mm             | 960x400x800mm          |
| Number of Baskets           | 2                     | 6                      |

